

Pharmacy Drawer Options

OmniceLL offers a wide variety of drawer options for storing, controlling, and accurately tracking the use of medications and pharmacy supplies. Pharmacy drawers integrate into the Omnicell OmniRx® automated dispensing cabinet product family and Anesthesia Workstation™ to create customized systems for specific clinical needs.

Modular Drawer Configurations

OmniceLL pharmacy drawer modules can be configured in a variety of ways within Omnicell automated dispensing cabinets. The OmniRx cabinet can accommodate up to 13 drawers, while the OmniTT™ tabletop cabinet can accommodate three pharmacy drawers.

Each cell of the OmniRx one-cell, two-cell, and three-cell cabinets can fit a three-drawer pharmacy module, a nine-drawer pharmacy module, or an OmniDispenser™ module (with or without pharmacy drawers). The pharmacy drawer modules can be placed in the middle and/or lower zone of the cell, to provide maximum flexibility in configuring the cabinet to meet the specific needs of a nursing area.

OmniDispenser Single-Dose Dispensing System

The OmniDispenser stores medications in a separate, secured location in the cabinet and dispenses each dose individually, so nurses only have access to the specific quantities requested. The system automatically tracks every dose, eliminating time-consuming count-backs to confirm remaining inventories of controlled substances. This is a significant time saver for nurses that gives them more time for patient care.

Key Features

- Eliminates the need to confirm remaining inventories of controlled substances (i.e. no count-backs), greatly reducing discrepancies.
- Push-button technology helps ensure accurate restocks.
- Security facilitates improved compliance with The Joint Commission and other regulatory standards.

Guiding Lights

OmniceLL's unique guiding lights technology visually directs clinicians to the drawer and then to the appropriate bin within the drawer for the medication they had selected. This feature reduces the possibility of nurses taking the wrong item.

FlexBin™ Single-Dose Drawer

The FlexBin single-dose drawer provides a high level of security for medications. With this drawer, access to a medication is restricted to the item and quantity chosen by the clinician. No other doses are available.

There is no need to confirm the remaining quantity or conduct end-of-shift counts, which results in fewer discrepancies.




Key Features

- Guiding lights feature illuminates the drawer location and individual product bins
- Single-deep drawer maximizes inventory in a minimum amount of space
- Three bin sizes available to accommodate a variety of medication types
- Bins can be utilized as single-dose or multi-dose
- Flexible drawer configurations: 96 small bins, 48 medium bins, and 48 small/large combinations

Locking Lid Drawer

Locking bins, suitable for controlled substances, provide a high level of security by restricting access to one preselected medication at a time. An audio alarm sounds if a clinician attempts to access a non-preselected compartment and notification is recorded and sent to pharmacy. The tamper-proof metal lids of each bin further ensure security.

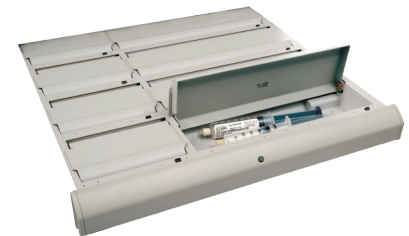

Key Features

- Guiding lights illuminate the drawer location and individual product bins
- Single-deep and double-deep drawers
- Four bin sizes: 6 bins, 8 bins, 12 bins, and 24 bins
- Tamper-proof, locked metal lids

Sensing Lid Drawer

Sensing bins, with see-through lids, provide a medium level of control. Special sensors on the bin lids detect bin entry and provide audible and visual feedback to clinicians attempting to open bins for items that have not been preselected. Any unauthorized entry attempts are recorded in the system, thus providing deterrence to inappropriate activity.



Key Features

- Guiding lights feature illuminates the drawer location and individual product bins
- Single-deep and double-deep drawers
- Three bin sizes
- Flexible drawer configuration: 6, 12, or 24 bins
- Clear plastic, entry-detecting, sensing lids
- Sensors provide audible and visual feedback
- Capacity up to 24 unique line items

Lighted Matrix Drawer

The Lighted Matrix Drawer provides a lower level of security but maximum flexibility. Lights frame the bin and guide the user to the preselected medication, reducing the possibility of errors.

Key Features

- Guiding lights illuminate drawer location and individual bins
- Flexible configuration with adjustable bin dividers
- Single-deep drawer
- Capacity up to 24 unique line items



Matrix Drawer

The Matrix Drawer is a high-capacity, low security drawer that is suitable for large quantities of non-controlled medications. The Double-Deep Matrix Drawer accommodates large quantities of medication or larger items such as large vials and premixed IV bags.

Key Features

- Adjustable bin dividers with capacity up to 96 unique items in a single-deep matrix drawer
- 4 non-adjustable bins in the double-deep drawer



Omnicell Medication/Supply Drawer

The Omnicell Medication/Supply Drawer provides a low level of security in width-configurable open bins. Placing this drawer within a pharmacy zone created in the cabinet can further control access to pharmacy items.



Key Features

- Guiding lights illuminate drawer location and individual product bins
- Easy-to-use green reorder buttons
- Single-deep drawer
- Adjustable bin dividers
- Capacity up to 24 unique line items

Nine-Drawer Pharmacy Module

The flexible design of the nine-drawer pharmacy module makes it the best solution for OmniRx one-cell, two-cell, and three-cell cabinets located in medication-intensive locations such as medical/surgical areas, emergency departments, and oncology areas, as well as transitional care and special care units.

Key Features:

- Perfect for areas where controlled substances are used
- Can be expanded for future needs: pay as you go for additional drawers
- Module accommodates many types of medications, including oral solids and liquids, syringes, ampules, vials, IV solutions, and PCA medications
- Can be configured with all drawer types to provide a complete solution for combined medication and supply management



Three-Drawer Pharmacy Module

The flexible, compact design of the three-drawer pharmacy module makes it ideal for locations that have a limited number of medications, such as the operating room, cath lab, and radiology, as well as urgent care and medical clinics.

- Accommodates many types of medications, including oral solids and liquids, syringes, ampules, vials, and bulk stock
- Can be configured with all Omnicell pharmacy drawer types